

**Better Vision  
Better Results**

## Firefly DE551

### Wireless Veterinary Video Otoscope



#### WIRELESS

With a dedicated wireless link and a 3 hour Lithium battery, DE551 provides complete flexibility and unprecedented precision in observing the ear canal and capturing high precision images / videos.

#### HIGH MAGNIFICATION

Provides up to 50x native optical magnification and up to 150x digital magnification

#### MULTI-LAYER GLASS LENSES

Equipped with multi-layer glass lenses rather than plastic alternatives, thus delivering superb focal depth and crystal clear image quality.

#### ULTRA-BRIGHT LEDs

Highest quality ultra-bright LEDs are deployed with variable brightness control

- **Performs highly accurate and clear otoscopy observations**
- **Also ideal for dental, dermal, nasal and rectal examinations**
- **Wirelessly captures snapshots & videos**
- **Magnifies objects up to 50x (optical) and up to 150x (digital)**

Firefly DE551 is the industry's leading Wireless high precision digital video otoscope with image and video capture capabilities. It provides unprecedented accuracy in observing the ear canal. It's also ideal for dental, dermal, nasal and rectal exams as well. This is a powerful tool which can be used for client education, veterinary training as well as electronic patient records. Unlike older electronic otoscopes, DE551 delivers these breakthrough capabilities - and much more - at a highly economical cost.




# Better Vision Better Results



## DE551

### Technical Details

Magnification	Native Optical: 15x – 50x Digital: 15x – 150x
Sensor Resolution	720 x 480
Lens Assembly	Dual Lenses 3-Layer Glass, 650nm cutoff
Video	Format: YUY2, MJPG Frame rate: 30 FPS
Video/Image Properties	Color: Hue, Saturation Exposure: Brightness, Contrast Image: Sharpness, Gamma
Image / Video Files	JPG / BMP / AVI
Battery	Built-in 850Mah
Illumination	4 Ultra-Bright LEDs Fully adjustable brightness
Wireless Link	4 Switchable Channels
Dimensions	13cm x 3.6cm x 4cm
Accessories	AC Battery Charger 3/4/5/7mm specula included
Interface	USB 2.0
Software	Scalable Window, Zoom, Freeze, Resolution, Rotate, Flip Region of Interest (ROI)
Operating Systems	Windows 7, 8, 10 and Mac OS-X 10.4 or higher
Warranty	1 Year Limited Hardware Warranty

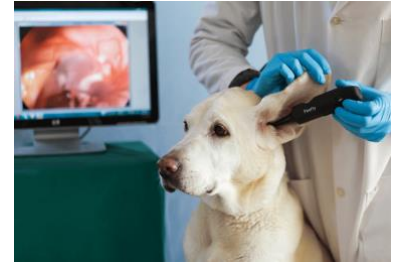

### SOPHISTICATED SOFTWARE

DE551 is controlled directly with the FireflyPro professional image processing software which is bundled with the product. This state-of-the-art software enables users to capture, store, recall, view, manipulate and measure images and videos in real time. Its intuitive interface empowers users to get right to work anytime and anywhere.



### Optional XL Specula Set (ES108)

This optional set includes six metal specula of various lengths including one with a working channel.

### Brought to you by:

### ADVANCED IMAGING TECHNOLOGY

Managing incoming light is a tremendous challenge in microscope design.

Reflections and lighting leaks can easily pollute the sensor, resulting in overexposure and washed-out images. Firefly effectively reduces overexposure by employing a proprietary lens assembly design and highly efficient image processing algorithms.

This highly efficient overexposure management technology enables Firefly to use ultra-bright LEDs to their fullest effect. Objects under observation are brightly illuminated without undesirable side effects.

Sophisticated image processing algorithms are ineffective without strong depth perception and a clear image. Firefly DE551 is equipped with multi-layer glass lenses rather than plastic alternatives, thus delivering superb focal depth and crystal clear images.

### INCLUDES WIRELESS RECEIVER



Specifications are subject to change in any matter and at any time without notice. The Firefly word and the Firefly logo are trademarks or trade names of Firefly Global and its affiliates in the United States and/or other countries. All rights in such names, marks or logos are reserved by Firefly Global.

Copyright © 2015 by Firefly Global, All rights reserved. Rev 1.2

[www.fireflyglobal.com](http://www.fireflyglobal.com)

**Better Vision  
Better Results**

A list of accessories compatible with the Firefly Veterinary Video Otoscopes.

## Firefly Video Otoscope Accessories

### ES108 XL Specula Kit



- 1pc 8mm diameter, 7cm long  
✓ Working Channel
- 1pc 7mm diameter, 7cm long  
✓ Working Channel
- 1pc 7mm diameter, 7cm long
- 1pc 6mm diameter, 6cm long
- 1pc 5mm diameter, 6cm long
- 1pc 4mm diameter, 5cm long

### ES140 Flexible Forceps



- Alligator forceps
- Flexible
- 13cm long
- Stainless steel
- Ideal for use with ES108

### ES160 Wireless TV Adapter



**Includes:**

- TV Output Cable
- AC Outlet Adapter
- Wireless Receiver

### ES101 Specula Set



**Includes three sizes:**

- 4mm (64 count)
- 5mm (64 count)
- 7mm (40 count)