

SISTEMA ELECTROQUIRÚRGICO ECO-800BI (PREMIUM)



SE PUEDE UTILIZAR PARA CIRUGÍAS
EN ANIMALES GRANDES.

ES ADECUADO PARA PROCEDIMIENTOS EN
ANIMALES PEQUEÑOS.

AYUDA A REDUCIR
EL SANGRADO INTRAOPERATORIO.

SE PUEDE APLICAR A PROCEDIMIENTOS
MÍNIMAMENTE INVASIVOS.

MORELIA, CDMX, GUADALAJARA, MÉRIDA, OAXACA, MONTERREY, HERMOSILLO, QUERÉTARO, PUEBLA, MIAMI, PANAMÁ, CHIRIQUÍ, COSTA RICA, COLOMBIA

☎ 443 233 0303 ✉ ventas@desego.com 🌐 desego.com 🇪🇸 🇫🇷 🇲🇽 🇮🇹 🇺🇸 **DESEGO**

FRECUENCIA DE SALIDA:

300/512KHZ

CORTE PURO:

300 W

CORTE MIXTO 1:

200 W

CORTE MIXTO 3:

80 W

COAGULACIÓN SUAVE:

120 W

COAGULACIÓN FORZADA:

120 W

COAGULACIÓN BIPOLAR 1:

80 W

MODOS DE CORTE



High Cut

Used for vaporization resection under water.



Pure Cut

Used for precise, non-damaging cutting of any tissue.



Blend Cut 1

Cutting mode with mild coagulation effect.



Blend Cut 2


Cutting mode with moderate coagulation effect.



Blend Cut 3

Cutting mode with stronger coagulation effect than Blend Cut 2.

MODOS DE COAGULACIÓN MONOPOLAR

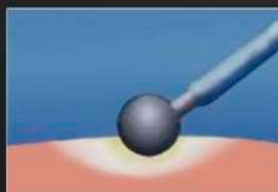

Soft Coagulation

Contact coagulation for hemostasis, promoting better postoperative wound healing.



Forced Coagulation

Rapid contact hemostasis and coagulation, also suitable for tissue cutting with coagulation.



Spray Coagulation

Suitable for hemostasis over large bleeding areas.